


ISI Surgical Instrument Cleaners and Lubricants, Cleaning Brushes and Brush Organizers

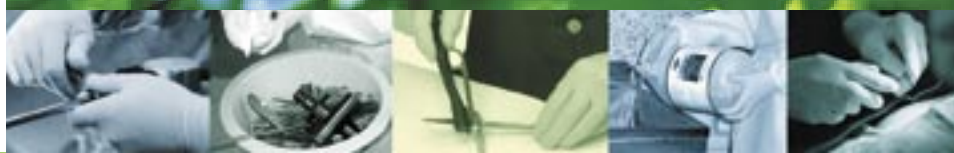

A Division of Instrument Specialists, Inc.

2007

1-800-537-1945 www.isisurgery.com
ISI for Ideas, Service, Integrity

Enzymatic Detergents
Detergents
Rust & Stain Removers
Lubricants
Instrument Brushes
Micro Brushes
ISI Brush Set Organizers

ISI Instrument Detergents, Cleaners, and Lubricants

Get the most from your surgical instruments with ISI instrument care solutions – FREE samples available



Back in the late 70's, our company's founder, Joseph R. Kabanek, Ph.D., noticed that many of the instruments sent for repair were not broken, but simply needed a thorough cleaning. His observation led to the development of our instrument care solutions – enzymatic detergents, detergents, rust and stain removers, and lubricants.

Today ISI continues to develop and perfect instrument care solutions to help keep delicate surgical and medical instruments functioning properly in hospitals across the country and around the world.

For your FREE samples, simply call our toll-free number or visit us at www.isisurgery.com

"...not broken, but simply needed a thorough cleaning..."



Free Samples 1-800-537-1945 www.isisurgery.com

ISI Quick Reference Guide

Effective, biodegradable and non-hazardous • Safe for plastic, rubber and metal

ISI Product	Ready to use	Concentrate [†]	For ultrasonic cleaner	For washer sterilizer	Manual cleaning	Soak
Enzymatic Detergents						
Hemaway Enzymatic Original		●	●	●	●	●
Hemaway Enzymatic Tri-Max		●	●	●	●	●
Hemaway Foam	●				●	●
EndoTerg, for flexible scopes		●		●	●	●
Detergents						
MPD, multi-purpose		●	●		●	●
Suds-Rite, low foam		●	●	●		
DuoKlenz, detergent, lubricant*		●	●	●		
Rust & Stain Removers						
SSK-80, stainless steel cleaner		●	●	●		
CSK-81, cool soak		●			●	●
Lubricants						
Instra-Glide Soak	●					●
Instra-Glide Jet		●		●		
DuoKlenz, lubricant, detergent*		●	●	●		

* Dual Purpose Cleaner † Concentrated formula makes 128 gallons of solution

■ Enzymatic Detergents

Effective, biodegradable and non-hazardous

Hemaway Original

Enzymatic strength removes organic materials

Hemaway Original removes dried blood, tissue, mucous, and other body fluids. Hard to reach places like box joints and tubing can't hide from the cleaning power of Hemaway Original. Reduces scrubbing time and prolongs the life of your instruments, Hemaway Original is formulated for daily use.

- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Removes organic materials, even dried-on blood
- ▶ Safe for plastic, rubber, metal

Hemaway Tri-Max

Hemaway Tri-Max contains the three strongest enzymes available for maximum cleaning

Similar to Hemaway Original, but formulated with the strongest combination of enzymes available. Concentrated Hemaway Tri-Max greatly reduces the need for scrubbing and prolongs the life of your instruments. Hemaway Tri-Max is strong enough to remove dried blood, tissue, mucous, and other body fluids, yet gentle enough to be used everyday.

- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Removes organic materials, even dried-on blood
- ▶ Maximum combination of enzymes available
- ▶ Safe for plastic, rubber, metal

Hemaway Enzymatic Foam

The ready-to-use, pre-cleaning foam for surgical instruments

Hemaway Enzymatic Foam is ideal for use in the O.R. to pre-clean surgical instruments and prevent corrosion. The foam contains active enzymes that go to work immediately to break down organic materials.

- ▶ Ready to use; spray directly on surgical instruments
- ▶ Removes organic materials safely
- ▶ Safe for plastic, rubber, metal

EndoTerg **NEW**

The low-foaming enzymatic detergent for flexible endoscopes

EndoTerg breaks down dried blood, tissue, mucous, and other bio-burdens on flexible endoscopes without leaving a soapy residue. Designed specifically for flexible endoscopes, EndoTerg provides thorough but gentle cleaning.

- ▶ Bacteriostatic and free rinsing
- ▶ Removes organic materials safely
- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Safe for plastic, rubber, metal

■ Detergents

Effective, biodegradable and non-hazardous

Duo-Klenz **NEW**

Removes blood and organic materials while lubricating surgical instruments

Highly concentrated, low-foaming Duo-Klenz removes blood and organic materials while safely lubricating surgical instruments – all in one easy step! Duo-Klenz is formulated for use in the wash cycle of washer sterilizers, cart washers, bottle washers, and other pressurized water units.

- ▶ Detergent combined with a non-silicone lubricant
- ▶ Removes organic materials, even dried-on blood
- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Safe for plastic, rubber, metal

MPD

The multi-purpose detergent

MPD is the highly concentrated multi-purpose detergent formulated to clean surgical instruments manually or with an ultrasonic cleaner. MPD is tough enough to remove blood and organic material and safe enough for everyday cleaning.

- ▶ For manual cleaning and use in ultrasonics
- ▶ Removes organic materials, even dried-on blood
- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Safe for plastic, rubber, metal

Suds-Rite

Low-foaming detergent concentrate for use in washer sterilizers and other pressurized water units

Use Suds-Rite in the wash cycle of washer sterilizers, cart washers, bottle washers, and other pressurized water units. Suds-Rite removes blood and organic materials from surgical instruments.

- ▶ The anti-foaming agent makes it perfect for use in the wash cycle of washer sterilizers and pressurized water units
- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Removes organic materials, even dried-on blood
- ▶ Safe for plastic, rubber, and metal



■ Rust & Stain Removers

Effective, biodegradable and non-hazardous

SSK-80

The stainless steel rust and stain remover for surgical instruments

SSK-80 is our best and most convenient instrument rust and stain remover for use in ultrasonic cleaners.

Used routinely, SSK-80 will extend the life of your instruments by removing and inhibiting rust, stains, and mineral deposits. The result: pivot points, box joints, and moving instrument parts function better. Simply add SSK-80 to your ultrasonic during the fill cycle. It's that easy!

- ▶ Concentrated; pour directly into the ultrasonic cleaner
- ▶ Removes rust, stains and mineral deposits
- ▶ Does not leave a residue
- ▶ Safe for plastic, rubber, metal

CSK-81

Concentrated rust and stain remover for soaking surgical instruments

Don't have an ultrasonic cleaner available? Or need to clean instruments not suitable for ultrasonic cleaners? Either way, CSK-81 instrument soaking solution removes rust, stains and mineral deposits leaving your surgical instruments cleaner and in better condition. Use routinely to inhibit new rust and extend the life of your instruments. Simply soak your surgical instruments and CSK-81 will do the rest!

- ▶ Concentrated formula makes 128 gallons of solution
- ▶ No heating required
- ▶ Removes rust, stains and mineral deposits
- ▶ Safe for plastic, rubber, metal



■ Lubricants

Effective, biodegradable and non-hazardous

Instra-Glide Soak

Non-silicone lubricant soak for surgical instruments

Instra-Glide Soak is the non-silicone lubricant that protects your surgical instruments. Instra-Glide Soak is water-soluble and provides a protective film without getting sticky or gummy. Keep your surgical instruments functioning at their optimum with this non-silicone lubricant.

- ▶ Ready to use; soak instruments prior to sterilization
- ▶ Water-soluble for easy rinsing
- ▶ Safe for plastic, rubber, metal

Instra-Glide Jet **NEW**

Non-silicone lubricant

Instra-Glide Jet is the non-silicone lubricant designed to keep surgical instruments functioning properly. Water-soluble Instra-Glide Jet provides a protective film without getting sticky or gummy. Helps prevent rust, stains and corrosion.

- ▶ Concentrated – makes 128 gallons of solution
- ▶ Water-soluble for easy rinsing
- ▶ Preservatives added to prolong shelf life
- ▶ Use in last cycle of washer to lubricate set

Duo-Klenz **NEW**

Removes blood and organic materials while lubricating surgical instruments

Highly concentrated, low-foaming Duo-Klenz removes blood and organic materials while safely lubricating surgical instruments – all in one easy step! Duo-Klenz is formulated for use in the wash cycle of washer sterilizers, cart washers, bottle washers, and other pressurized water units.

- ▶ Detergent combined with a non-silicone lubricant
- ▶ Removes organic materials, even dried-on blood
- ▶ Concentrated formula makes 128 gallons of solution
- ▶ Safe for plastic, rubber, metal



**Request Your
Free Samples**
1-800-537-1945
www.isisurgery.com

Brush Sets & Organizers

ISI Brush Nest, ISI Micro Nest wall-mount organizers

Get organized for effective instrument cleaning

The ISI Brush Nest and Micro Nest wall organizers come with a complete set of cleaning brushes in a variety of sizes and diameters. Both the Brush Nest and Micro Nest are designed to help technicians select the best brush for the cleaning task.

Individual packages of replacement brushes are available for both the ISI Brush Nest and ISI Micro Nest.



ISI Brush Nest

Made of tough durable plastic with a convenient self-draining shelf, the ISI Brush Nest comes with 18 specialized brushes in different lengths, diameters and compositions.



ISI Micro Nest

Made of durable plastic, the ISI Micro Nest organizes a complete set of micro brushes. All micro brushes have stainless steel handles and nylon bristles, plus come in a variety of sizes and French diameters. These brushes are designed for hard-to-clean instruments such as cannulas.



Brushes

ISI brushes are durable and come in a variety of diameters and sizes



Instrument Cleaning Brushes

ISI offers durable instrument cleaning brushes in a variety of lengths, diameters and compositions. Sold in sets and individually.

Micro Brushes

ISI Micro Brushes are available in a variety of French diameters. Made with stainless steel handles and nylon bristles, ISI Micro Brushes will last through countless cleanings. Sold in sets and individually.

ISI Free Offers

Call ISI at 1-800-537-1945 or visit www.isisurgery.com

Free samples of ISI instrument cleaners and lubricants

If you haven't tried our instrument cleaners and lubricants, here's your chance! Call or go online for your free samples.

Free Alumi-Hand samples

With the Alumi-Hand you trade the hazards of lead for a sterile, disposable, affordable surgical hand immobilizer. Call or go online for your free sample.

Free ISI surgical table evaluations

Don't just take our word that these are great surgical platforms – see for yourself! The Universal K, Advanced Universal K, Self-Supporting Arm Table, Ultra-Light, and Free-Standing Arm Table are available for a 30 day, no-charge evaluation in your O.R.

Two divisions, the same commitment

Since 1978, the professionals at ISI have worked together with pride. A unique dedication to craftsmanship and quality at the best possible price has earned ISI trust across the medical community.

ISI Medical & O.R. Products

- ▲ Specialty Surgical Tables
- ▲ Leadless Hand Immobilizers
- ▲ Upper Extremity Traction Devices
- ▲ Instrument Cleaning Solutions
- ▲ Instrument Lubricants
- ▲ Instrument Cleaning Brushes, Sets and Brush Organizers
- ▲ GSA's Federal Supply Service (FSS)
- ▲ TXMAS Contract Pricing
- ▲ Other GPO Pricing Available

ISI Medical Equipment Service & Sales

Service and Repairs:

- ▲ Flexible & Rigid Endoscopes
- ▲ Medical Camera Repair
- ▲ Monitor & Printer Repair
- ▲ Powered Hand Piece Repair
- ▲ Phaco Systems Repair
- ▲ Endoscopic Accessories
- ▲ GSA's Federal Supply Service (FSS)
- ▲ TXMAS Contract Pricing
- ▲ Other GPO Pricing Available

ISI, always at your service

We look forward to the opportunity to serve our clients as well as earn the respect of new customers. For more information visit www.isisurgery.com or speak with one of our customer representatives at 1-800-537-1945.

Instrument Specialists, Inc.

www.isisurgery.com ▲ Toll Free in the USA: (800) 537-1945

Tel: (830) 249-9535 ▲ Fax: (830) 249-9433

32390 IH-10 West, Boerne, Texas 78006-9214

