

# ISI Lexan-Top Surgical Tables – For Upper Extremity Surgeries

## Clearly a good choice for hand surgeries and upper extremity cases

ISI's Lexan-top tables have been around for years and the reason is clear – they're great surgical tables for upper extremity cases. So if you're looking for extraordinarily strong and dependable hand tables, look no more. The answer is clear: ISI's Lexan Tables.

### The Self-Supporting Arm Table



**The Self-Supporting Arm Table** is the ISI's "legless" table. Built-in mounting brackets attach to the side rail of your OR Table to provide the support needed for a table that allows easy access for surgeons, staff and equipment. An optional leg is available if needed.



The Self-Supporting Arm Table is designed for easy access for C-arm and staff.

### The Lexan Universal K Table

**The Lexan Universal K Table** has a sturdy, two-prong leg that creates the steadiest operative platform available. This surgical table attaches to the side rail of your OR Table and has the option of an Upper Extremity Traction System. An optional mounting piece is available for attaching to the end of your OR table for head & neck or foot & ankle surgeries.



The Lexan Universal K Table optional upper extremity traction system



**ISI Medical & O.R. Products**

A DIVISION OF INSTRUMENT SPECIALISTS, INC.

1-800-537-1945 [www.isisurgery.com](http://www.isisurgery.com)

## The Self-Supporting Arm Table



### Specifications

- STOCK NO.** AT-4000
- DIMENSIONS** 20" wide by 30" long (curved inset in table top, 16-1/2" wide)
- WEIGHT** 18 pounds
- MATERIALS** G.E. Lexan top, stainless steel mounting brackets

### OPTIONS

- Optional Single-Post Telescoping Leg
- 2" and 4" contoured foam pads



### Care Instructions

Wash table and pad with soap and water and/or disinfectant. Do not clean with abrasive cleanser.

**The Self-Supporting Arm Table is designed for easy access for C-arm and staff**



## The Lexan Universal K Table




### Specifications

- STOCK NO.** K3510
- DIMENSIONS** 22-1/4" wide by 36" long (curved inset in table top, 18-1/2" wide)
- WEIGHT** 20 pounds
- MATERIALS** G.E. Lexan top, Aluminum support struts, stainless steel tubing for leg and mounting brackets
- OPTIONS**
- 2" and 4" contoured foam pads
  - Traction package of finger traps and weights
  - Choice of mounting bracket

### Care Instructions


Wash table and pad with soap and water and/or disinfectant. All parts of the traction set may be sterilized, with the exception of the carrying case.

### Partial list of the various types of surgeries performed on the Lexan Universal K Table

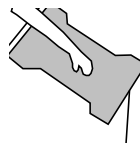


#### EXTREMITY SURGERIES

- Wrist Arthroscopies
- Ankle Arthroscopies
- Carpel Tunnel Release
- Fracture and Pin Replacement
- Amputations
- Grafts, Reconstructive and Plastic Surgeries

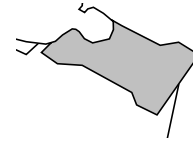


#### HEAD/NECK SURGERIES

- Carotid Angiogram & Endarterectomy
- Maxillofacial Surgeries
- Laryngostomy
- Grafts, Reconstructive and Plastic Surgeries

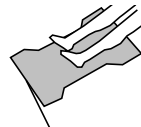

#### GENERAL SURGERIES

- AV Shunts
- Venous Cutdowns
- Tumor Removal
- Bones Grafts

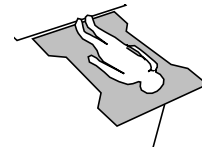

#### PEDIATRIC SURGERIES

- Orthopedic Surgeries, i.e. Hip Dysplasia
- Arthrograms and Angiograms
- Grafts, Reconstructive and Plastic Surgeries
- General Pediatric Surgeries



**ISI Medical & O.R. Products**

A DIVISION OF INSTRUMENT SPECIALISTS, INC.

1-800-537-1945 [www.isisurgery.com](http://www.isisurgery.com)